



max surgical

Distribuidor exclusivo para o Brasil

Cambria™



Cambria™

Anterior Cervical Interbody System

Cambria™ Anterior Cervical Interbody System




The Cambria Cervical Interbody System combines unique dual convex surface geometry, a wide variety of implant options, and a comprehensive set of instrumentation to limit uncertainty for anterior cervical interbody fusion.

Anatomic Solution

- Contoured surface geometry matches that of the endplates to maximize engagement.
- Tapered leading edge eases implant insertion.
- Large graft window to aid in fusion.

PEEK-OPTIMA®

The properties of PEEK-OPTIMA® thermoplastic are optimal for spinal implants. Its mechanical strength is close to that of cortical bone. PEEK is a biocompatible, radiolucent polymer which allows clear assessment of bony fusion.

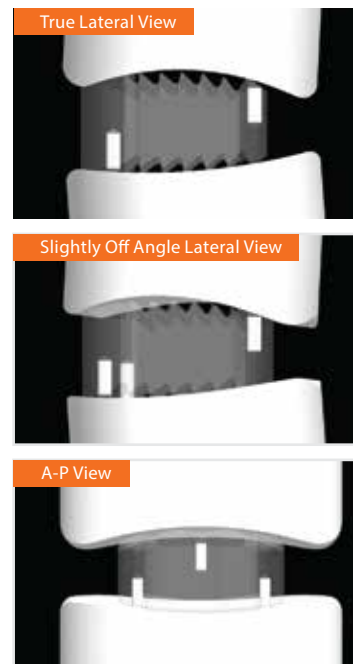

Tantalum markers assist in the visualization of the Interbody position in the implant site.

Cambria Design Advantages

- Three footprint options: 13mm, 15mm, 17mm widths
- Implant Trials available with and without depth limiting stops.
- Unique tooth pattern reduces risk of migration.


6.5 degrees built in lordosis with contoured superior and inferior surfaces

Central Graft window provides large area for fusion


Manufactured from PEEK-OPTIMA® for superior biocompatibility and modulus close to that of natural bone



Trials with depth limiting stops

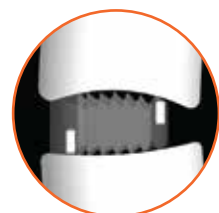


Dual Convex Surface

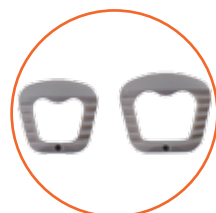

Multiple footprint options



Retentive Inserter enables accurate insertion

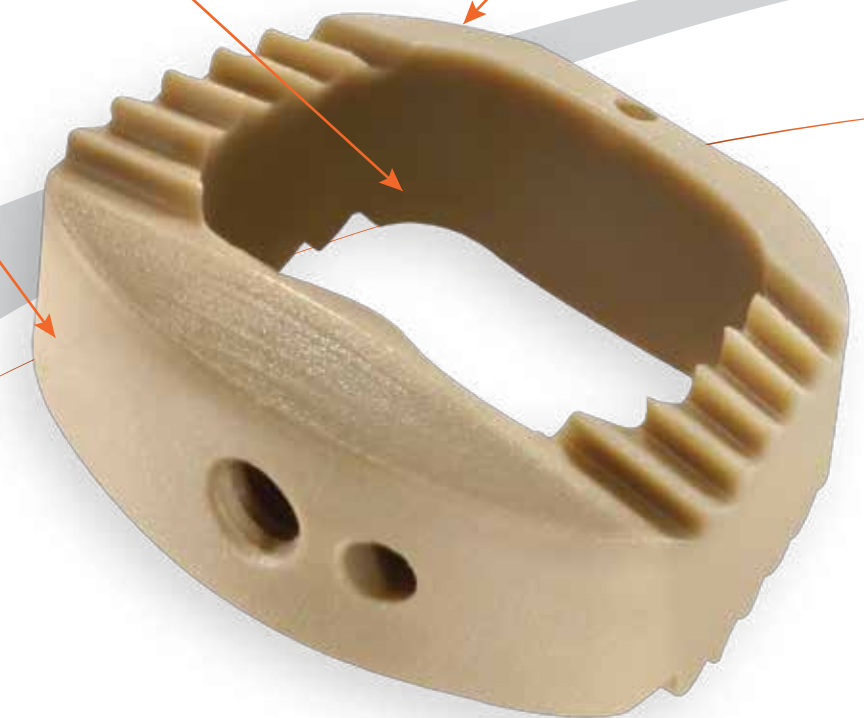

System Specifications

Cambria Implant Sizes (6.5° and Flat*) (Depth x Width x Height)

12 x 13 x 5mm	13 x 15 x 5mm	13.5 x 17 x 5mm
12 x 13 x 6mm	13 x 15 x 6mm	13.5 x 17 x 6mm
12 x 13 x 7mm	13 x 15 x 7mm	13.5 x 17 x 7mm
12 x 13 x 8mm	13 x 15 x 8mm	13.5 x 17 x 8mm
12 x 13 x 9mm	13 x 15 x 9mm	13.5 x 17 x 9mm
12 x 13 x 10mm	13 x 15 x 10mm	13.5 x 17 x 10mm
12 x 13 x 11mm	13 x 15 x 11mm	13.5 x 17 x 11mm
12 x 13 x 12mm	13 x 15 x 12mm	13.5 x 17 x 12mm

*Flat implants are available upon request.



Para mais informações ou para fazer um pedido, entre em contato:

max surgical

Distribuidor exclusivo para o Brasil

T: +55 (19) 3522-2550
+55 (19) 3524-6964

Avenida 4, 1317 - Jardim Claret
CEP 13503-230 - Rio Claro, SP, Brasil

www.maxsurgical.com.br

Manufacturer:



SeaSpine, Inc.
5770 Armada Drive
Carlsbad, CA 92008 USA
Phone 866.942.8698 USA | Fax 760.525.0337



mdi Europa GmbH
Langenhagener Str.71
30855 Hannover-Langenhagen Germany
Phone: +49 511 39089530 Fax: +49 511 39089539
info@mdi-europa.com www.mdi-europa.com

Warning: Applicable laws restrict these products to sale by or on the order of a physician.

Cambria and the SeaSpine logo are trademarks of SeaSpine Orthopedics Corporation or its subsidiaries in the United States and/or other countries. ©2016 SeaSpine Orthopedics Corporation. All rights reserved.